

Kowa



Technology for Life Science

KOWA APPLANATION
TONOMETER


HA-2



Convenient Item For Daily Clinical Practice



● Standard set of HA-2

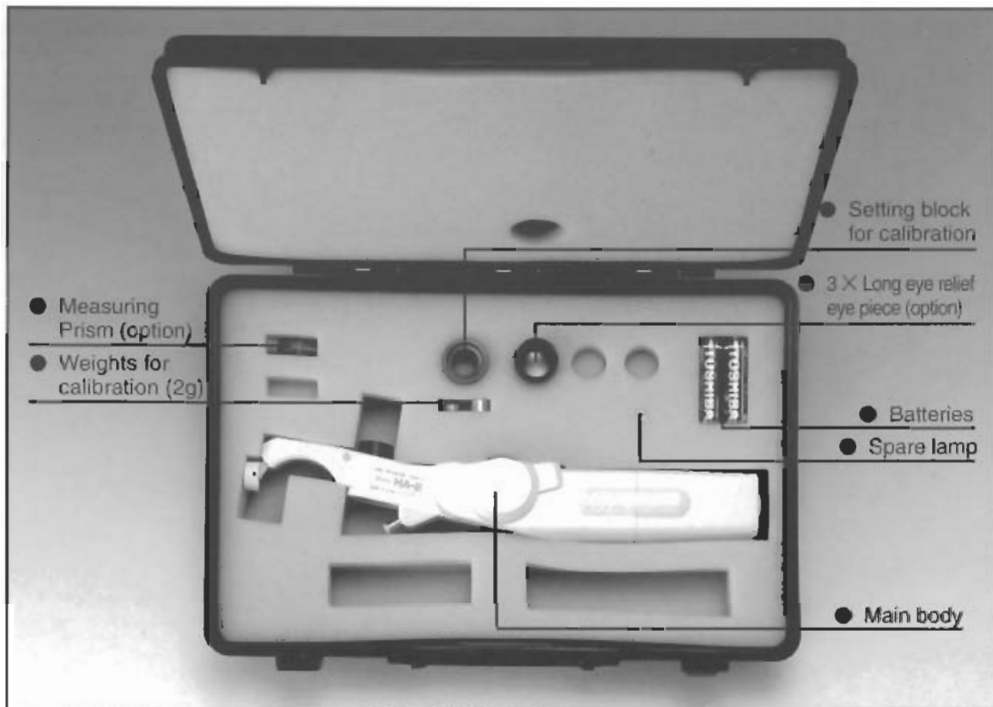

3 X Long eye relief eye piece (option)

● Main Features

- ◆ Compact and light-weight, well balanced and easy to handle by single hand.
- ◆ Unique one-spring mechanism assures correct reading regardless of patient posture.
- ◆ Built-in illumination permits you to use the instrument regardless of the illuminating condition of consultation room.
- ◆ Direct reading type precision dial requires no conversion table to change the readings.
- ◆ Accessories such as adjustable head rest and interchangeable eyepiece provide most suitable conditions for measurement.

● Specifications

Type	Hand-held applanation tonometer
Applanation prism	Doubling prism interchangeable (option)
Scale range	0-60 mm Hg (1 mm Hg division)
Illumination	
a) ON/OFF	Interlocked with thumb wheel
b) Bulb	Lens bulb 2.5 V, 0.25 A interchangeable
c) Battery used	1.5 V, 2 pcs Interchangeable
Eyepiece Magnification	6 X
Head rest	Movable
Length	290 mm
Weight	240g with batteries

* All specifications and external design are subject to modification of change without prior notice.



Notified Body No.0197



Kowa Company, Ltd. has been certified for quality system in accordance with ISO9001, EN46001.


Kowa Kowa Company, Ltd.

Life Science Division
No.4-14, 3-chome, Nishinbashi-Hondcho, Chuo-ku, Tokyo 103-8433, Japan
Phone : 81-3-3279-7311
Facsimile : 81-3-5255-7516

Kowa Europe GmbH *Kowa Optimed, Inc.*

Limmenmanntstrasse 65A
40210 Düsseldorf F.R.G. Germany
Phone : 49-211-1793540
Facsimile : 49-211-161952
URL : <http://kowa-europe.com/>

20001 S. Vermont Ave. Torrance,
CA 90502, U.S.A.
Phone : 1-310-327-1913
Facsimile : 1-310-327-4177
URL : <http://medical.kowa.com/>

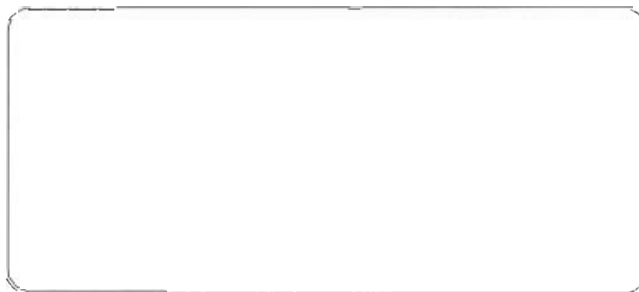

PRINTED WITH SOY INK This catalog was made with 100% recycled paper

Printed in Japan 0389500MA